




Hydro-Trach T Mk.II

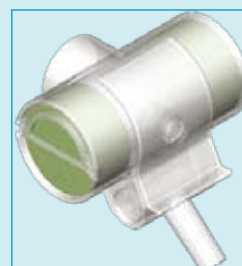

Hydro-Trach T Mk.II

Designed for use on tracheostomised patients


The Hydro-Trach T[®] Mk.II is a heat and moisture exchange device designed for use with spontaneously breathing patients in order to reduce loss of heat and moisture during respiration.

When a patient has a tracheostomy, the normal system of temperature and moisture maintenance is bypassed by the insertion of the tracheal tube. The possible loss of heat and moisture can lead to serious complications, notably damage to cilia and the mucous glands.

This in turn may result in retention of sputum and atelectasis, production of mucous plugs and potential tubal occlusion.

The Hydro-Trach T[®] Mk.II has a number of unique features which make it an ideal product for prolonged use with spontaneously breathing patients.

Features and benefits

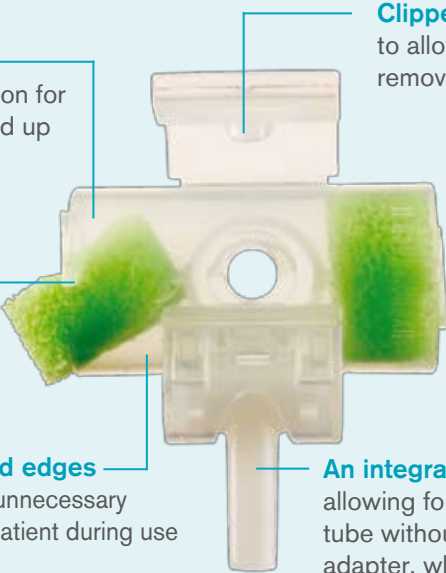

- Clear housing** for easy visual inspection for possible secretion build up
- Anti Occlusion mechanism** Safety mechanism allowing the HME element to partially dislodge in the event of total occlusion or vigorous cough, allowing ease of breathing
- Clipped suctioning port** to allow easy suctioning without removal of the device
- Standard 15mm fitting** fits securely onto standard tracheal tubes
- Small and lightweight** reducing the pull on the patient connection
- Smooth round edges** to prevent the unnecessary trauma to the patient during use
- An integral swivel oxygen connector** allowing for easy connection of the oxygen tube without the need of a separate oxygen adapter, which can be easily folded away when not in use

Technical Specification

Moisture return at: VT 500ml HME only	Resistance at: HME only		Compressible volume HME only	Weight HME only	Minimum tidal volume HME only
	30L/min	60L/min			
26mg H ₂ O/L	0.2cm H ₂ O	0.7cm H ₂ O	19ml	8g	50ml

Ordering information

Code: 1873000	Box qty. 25
Hydro-Trach T Mk.II HME	
Code: 1874000	Box qty. 40
Hydro-Trach T Mk.II + O ₂ tube	
Code: 1873000S	Box qty. 100
Hydro-Trach T Mk.II HME sterile	
Code: 1874000S	Box qty. 30
Hydro-Trach T Mk.II + O ₂ tube sterile	

IS12.2 • Issue 2 11.08