



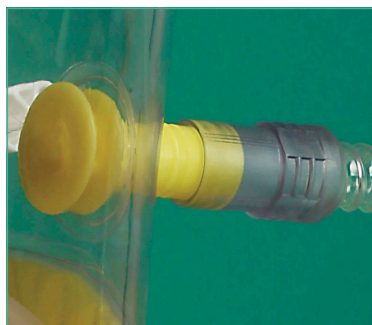

STARMed CaStar R

Hood for non-invasive mechanical ventilation (NIV)

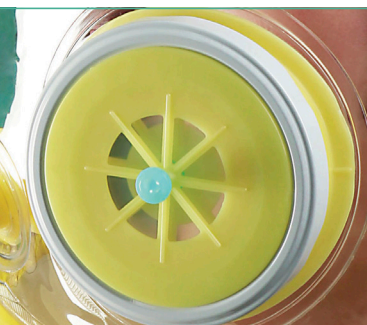

Critical Care ▪ Non-invasive ventilation ▪ **STARMed**



STARMed CaStar R hood for NIV

StarMed CaStar R

The CaStar R is a comfortable, versatile and lightweight patient interface for non-invasive mechanical ventilation (NIV) designed to offer an alternative to traditional methods of NIV delivery. The CaStar R is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in both semi-recumbent and supine positions providing an ideal option for NIV therapy.

Technical Features

- 22M connectors compliant to ISO standard
- Internal volume 11 litres - weight 380g
- Comfortable and atraumatic underarm fastening straps
- Built in inflatable cuff

- Two sealed access ports for probes or catheters of 3.5-7.0mm Ø
- Patient access port with screw cap. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

Bibliography:

- Mojoli F, et al. An optimized set-up for helmet noninvasive ventilation improves pressure support delivery and patient-ventilator interaction, Intensive Care Med. 2013 Jan;39(1):38-44
- Vargas F, et al. Helmet with specific settings versus facemask for noninvasive ventilation, Crit Care Med. 2009 Jun;37(6):1921-8
- Mojoli F, et al. Carbon dioxide rebreathing during non-invasive ventilation delivered by helmet: a bench study, Intensive Care Med. 2008 Aug;34(8):1454-60
- Antonelli M, et al. New treatment of acute hypoxemic respiratory failure: Noninvasive pressure support ventilation delivered by helmet - A pilot controlled trial, Crit Care Med. 2002 Mar;30(3):602-8

Features and benefits

Patient access port

With bi-directional anti-asphyxiation valve


Sealed access ports

For probes and catheters

Transparent

For ease of patient monitoring

Comfortable underarm fastening

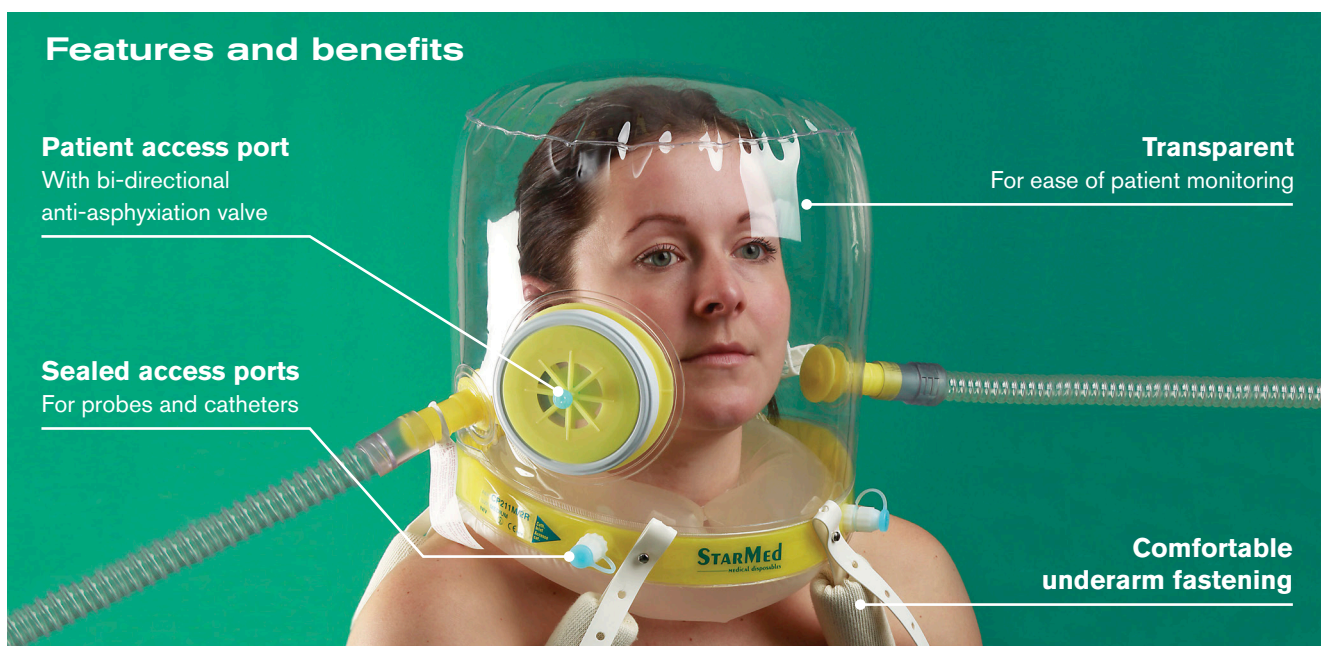

CaStar R hood for NIV with bi-directional anti-asphyxiation valve and patient access port

Custom packs are supplied on specific request of the customer

Code	Size	Neck circumference	Box Qty.
CP211XXL/2R	Extra extra large	≥50cm	5
CP211XL/2R	Extra large	45-52cm	5
CP211L/2R	Large	40-47cm	5
CP211M/2R	Medium	34-41cm	5
CP211S/2R	Small	27-34cm	5
CP211XS/2R	Extra small	17-27cm	5



www.intersurgical.com/info/starmed



CE certified product

Manufactured by: Intersurgical S.p.A.

IS1.18_ROW • Issue 4 03.15



Manufactured by: Intersurgical S.p.A. – Via T. Morandi, 12 – 41037 – Mirandola MO – Italy

T: +39 0535 610131 F: +39 0535 610310 info@intersurgical.it www.intersurgical.com

UK • France • Germany • Spain • Portugal • Italy
Netherlands • Sweden • Lithuania • Russia • South Africa
Japan • Philippines • USA • Colombia • Australia



Please think before you print
Save energy and paper.
If you must print this document
please print it double sided

Intersurgical Ltd is certified to ISO 14001, ISO 9001 and ISO 13485