




# STARMed CaStar R

Hood for non-invasive mechanical ventilation (NIV)

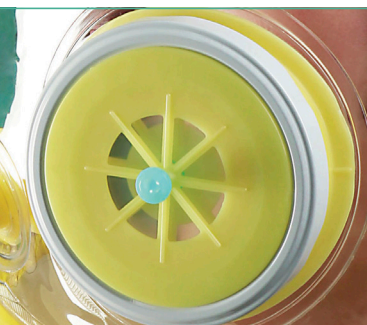

Critical Care ▪ Non-invasive ventilation ▪ **STARMed**



# STARMed CaStar R hood for NIV

## StarMed CaStar R

The CaStar R is a comfortable, versatile and lightweight patient interface for non-invasive mechanical ventilation (NIV) designed to offer an alternative to traditional methods of NIV delivery. The CaStar R is well tolerated for extended periods and can reduce the risks associated with both endotracheal intubation and face mask use. It can be used in both semi-recumbent and supine positions providing an ideal option for NIV therapy.


### Technical Features

- 22M connectors compliant to ISO standard
- Internal volume 11 litres - weight 380g
- Comfortable and atraumatic underarm fastening straps
- Built in inflatable cuff

- Two sealed access ports for probes or catheters of 3.5-7.0mm Ø
- Patient access port with screw cap. Fitted with bi-directional anti-asphyxiation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated

#### Bibliography:

- Mojoli F, et al. An optimized set-up for helmet noninvasive ventilation improves pressure support delivery and patient-ventilator interaction, Intensive Care Med. 2013 Jan;39(1):38-44
- Vargas F, et al. Helmet with specific settings versus facemask for noninvasive ventilation, Crit Care Med. 2009 Jun;37(6):1921-8
- Mojoli F, et al. Carbon dioxide rebreathing during non-invasive ventilation delivered by helmet: a bench study, Intensive Care Med. 2008 Aug;34(8):1454-60
- Antonelli M, et al. New treatment of acute hypoxemic respiratory failure: Noninvasive pressure support ventilation delivered by helmet - A pilot controlled trial, Crit Care Med. 2002 Mar;30(3):602-8

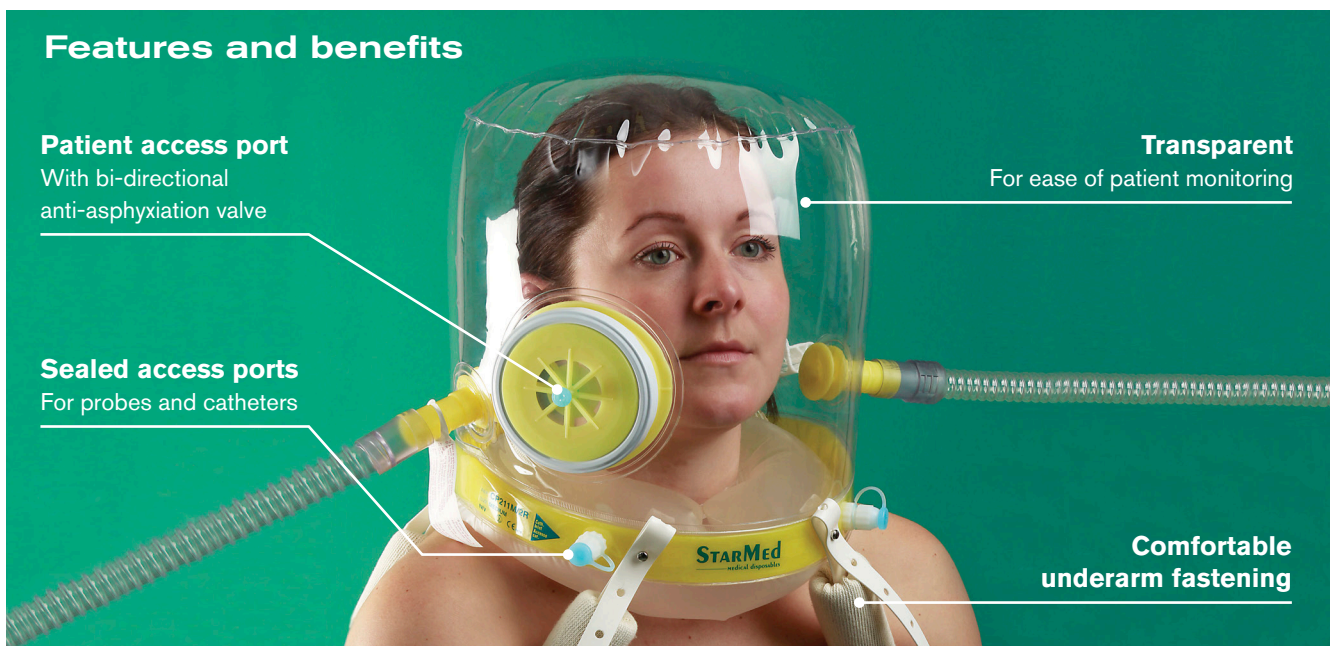

CaStar R hood for NIV with bi-directional anti-asphyxiation valve and patient access port, inflation bulb, ear plugs and measuring tape

Code	Size	Neck circumference	Box Qty.
CP211XXL/2RW	Extra extra large	≥50cm	5
CP211XL/2RW	Extra large	45-52cm	5
CP211L/2RW	Large	40-47cm	5
CP211M/2RW	Medium	34-41cm	5
CP211S/2RW	Small	27-34cm	5
CP211XS/2RW	Extra small	17-27cm	5



[www.intersurgical.co.uk/info/starmed](http://www.intersurgical.co.uk/info/starmed)



CE certified product

Manufactured by: Intersurgical S.p.A.

IS1.18\_UK • Issue 4 09.14



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
T: +44 (0)118 9656 300 F: +44 (0)118 9656 356 info@intersurgical.com www.intersurgical.com

UK ▪ France ▪ Germany ▪ Spain ▪ Portugal ▪ Italy ▪ Netherlands  
Sweden ▪ Lithuania ▪ Russia ▪ South Africa ▪ Japan ▪ Philippines  
USA ▪ Australia



**Please think before you print**  
Save energy and paper.  
If you must print this document  
please print it double sided

Intersurgical Ltd is certified to ISO 14001, ISO 9001 and ISO 13485