

STARMed CaStar Up

Openable hood for CPAP therapy



STARMed CaStar Up – openable CPAP hood

StarMed CaStar Up

Openable CaStar Up is a device recommended for the treatment of respiratory failure:


- Cardiogenic pulmonary oedema
- Non cardiogenic pulmonary oedema (ALI-ARDS)
- Thoracic trauma
- Atelectasis
- Post operative hypoxemia
- Asthma
- COPD exacerbation
- Pre and post extubation respiratory support

Advantages

- Allows rapid access to the patient's face during nursing procedures or in case of emergency
- One operator can easily remove the upper part of the hood and have full access to the patient
- Probes and catheters can remain in place, on the patient, connected to the lower part of the hood
- Well tolerated for extended periods
- Prevents the risk of pressure sores from face masks

Technical Features

- 22M/22F connectors compliant to ISO standard
- Fitted with integrated pressure gauge to check the presence of pressure within the hood
- Internal volume 15 litres – weight 690g
- Comfortable underarm fastening straps
- Other fastening accessories are available if required
- Two new larger access ports on the collar for probes or catheters of 3.5-7.0mm Ø
- Patient access port with screw cap. Fitted with bidirectional anti-suffocation valve which opens automatically in the event of pressure failure
- Airtight sealing of the system is achieved by an elastic membrane making contact with the upper part of the torso, once the hood is inflated
- Additional models available with a built-in inflatable neck cushion or without patient access port

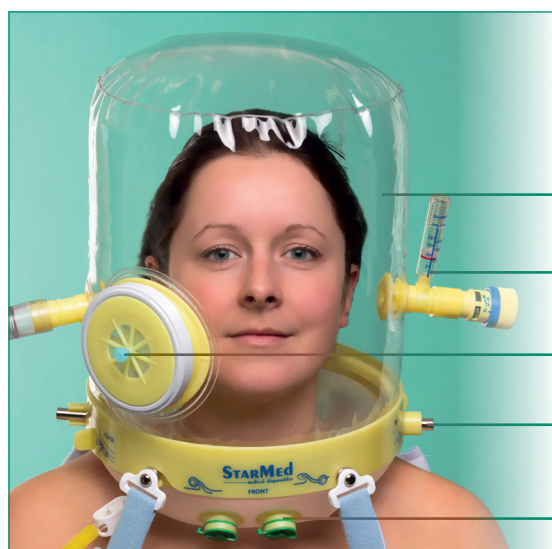

Features and benefits



CE certified product

Transparent

For ease of patient monitoring

Integrated pressure manometer

Patient access port with anti-suffocation valve

Easy removal for rapid patient access

New large dimension accesses on the collar

Ideal for single and multiple lumen probes or catheters

Castar Up CPAP hood with anti-suffocation valve and patient access port

Code	Size	Neck Circumference	Box qty.
CA011XXL/2	Extra Extra Large	≥50cm	5
CA011XL/2	Extra Large	45-52cm	5
CA011L/2	Large	40-47cm	5
CA011M/2	Medium	34-41cm	5
CA011S/2	Small	27-34cm	5
CA011XS/2	Extra Small	17-27cm	5

Single packed. Minimum purchasing package: 5 units. Custom-packs are supplied on specific request of the customer.



Manufactured by:
Intersurgical S.p.A. - Via T. Morandi, 12
41037 Mirandola MO - Italy
T. +39 0535 610131 - F. +39 0535 610310

UK
info@intersurgical.co.uk
Italy
info@intersurgical.it

España
info@intersurgical-es.com
Portugal
info@intersurgical.pt

Lietuva
info@intersurgical.lt
Россия
info@intersurgical.ru

France
info@intersurgical.fr
South Africa
info@intersurgical.co.za

Japan
info@intersurgical.co.jp
Nederland
info@intersurgical.nl

IS1.31_ROW • Issue 2_01.14

