

Nasal Cannulae

A range of cannulae for adults, paediatrics and neonates



Nasal Cannulae

Intersurgical is pleased to offer a complete range of nasal cannulae. This comprehensive range is entirely latex-free and includes designs suitable for use with adults, paediatrics and neonates.

Choice of options


Product choices are available to suit neonatal, paediatric and adult patients. The following prong shapes are also available:

- Curved
- Flared
- Curved/Flared
- Straight


All oxygen tubing in the range is manufactured in a 'multi-channel' design, which permits oxygen flow even if the tubing is 'kinked'.

Now available with ear guards


Ear guards are an additional comfort benefit that come pre-attached to the cannulae designed for patient satisfaction.

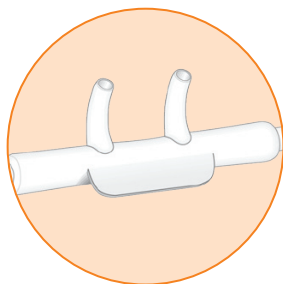


Curved prongs:
Improved anatomical fit

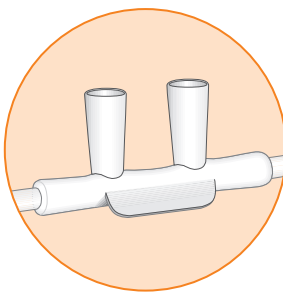

Flared prongs:
Slows down the flow of oxygen

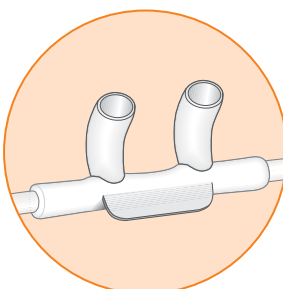

Curved / Flared prongs:
Combines the benefits of both the above designs

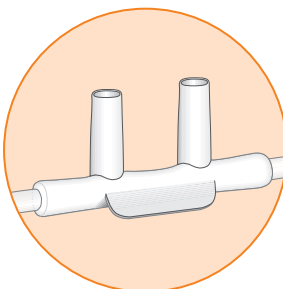

Straight prongs

Features and Benefits

Multi-Channel oxygen tubing

Allows oxygen to continue to flow if the tubing becomes kinked

Latex free

Suitable for patients with latex allergy

Over-ear design

Secure positioning
(Ear guards available for added comfort)

Non-slip toggle

Enables user to adjust for best possible fit

Clear material

Unobtrusive

Soft prongs

Minimise irritation

Location saddle


Assists in the optimum positioning of the prongs

Adult headset version


Ideal for domiciliary use




Ordering information

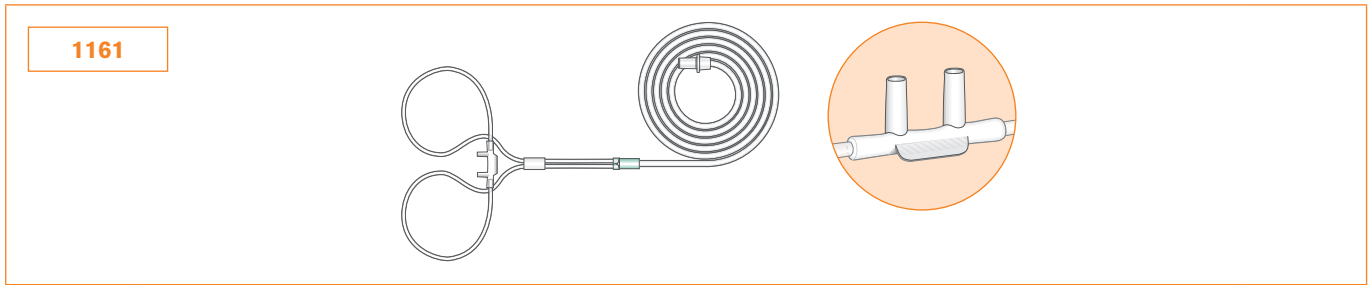


1161	Adult Straight Prong with tube, 1.8m length	50	
1162	Adult Straight Prong with tube, 5.0m length	20	
1169	Adult Straight Prong headset, 0.5m length	100	

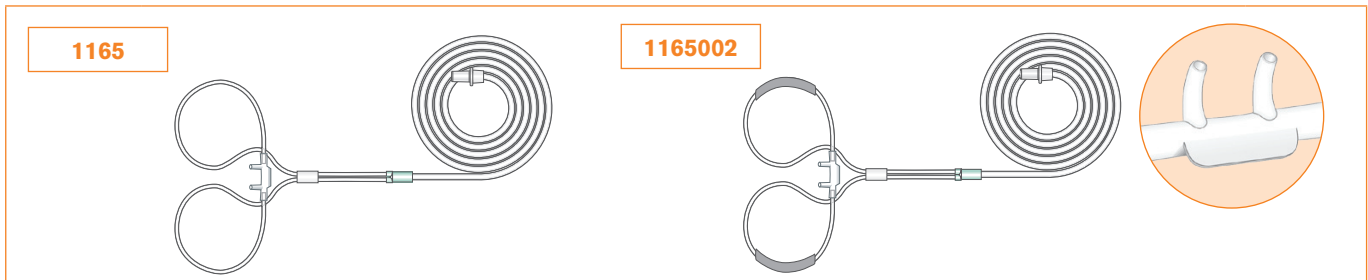

1165	Adult Curved Prong with tube, 1.8m length	50	
1165002	Adult Curved Prong with ear gaurds, 2.1 m length	50	
1163	Paediatric Curved Prong with tube, 1.8m length	50	
1164	Neonatal Curved Prong with tube, 1.8m length	50	
1165002	Adult Curved Prong with ear guards and tube, 1.8m length	50	

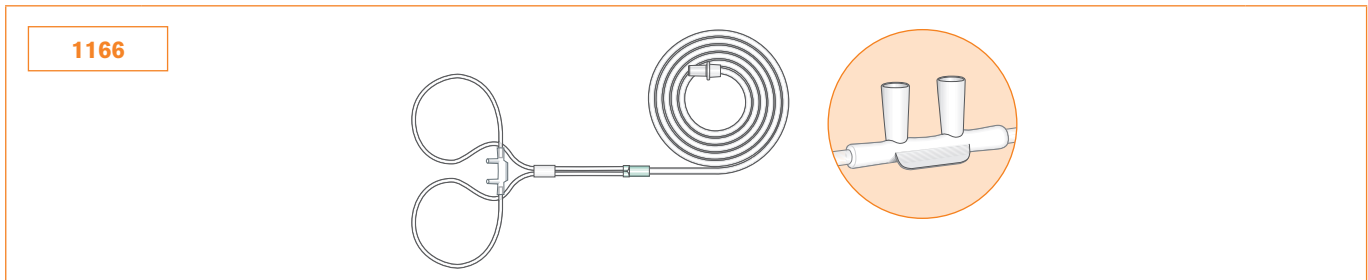

1166	Adult Flared Prong with tube, 1.8m length	50	
-------------	---	----	--

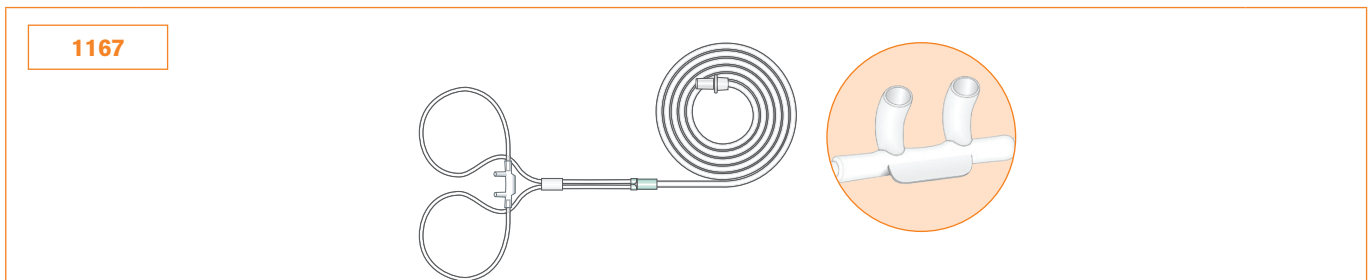

1167	Adult Curved / Flared Prong with tube, 1.8m length	50	
1168	Adult Curved / Flared Prong headset, 0.5m length	100	

IS10.8 • Issue 4 01.12